


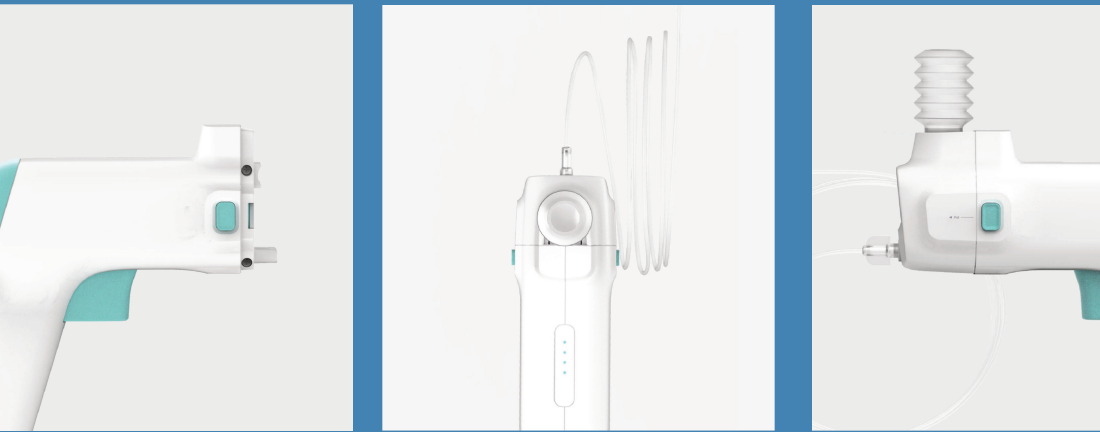

Customer Service

Official Website: www.lbbc.com

Phone: +86 18323413839

Email: services@cqlbbc.com

Whatsapp:



Yao




Meixi Wang

XunNing Portable Powder Spray Unit 2.0


A better Tool for endoscopic surgical hemostasis

What is XunNing Portable Powder Spray Unit 2.0 ?

An endoscopic spray delivery device specifically designed for use with Xunning Hemostatic Powder. Through a pistol-style air pump tube and delivery catheter, it enables rapid and precise spraying of Xunning Hemostatic Powder onto bleeding sites within the body via the endoscope channel, achieving efficient hemostasis.

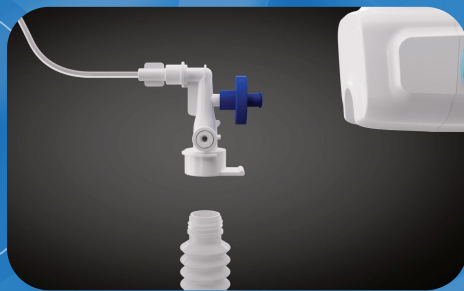


①Position the gas-powder mixer interface vertically downward, align it with the Xunning hemostatic powder bottle, after tightening, flip the mixer.

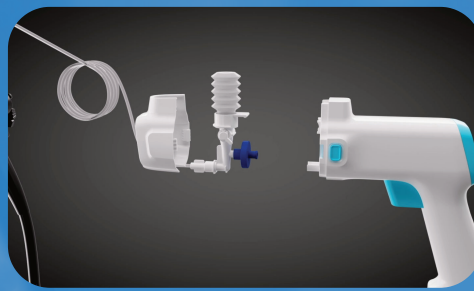


②Insert the mixer into the spray device's positioning connector, and make sure it is securely clamped.



③Turn on the power switch at the rear of the spraying device.

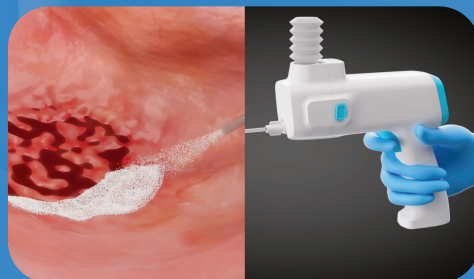

④Position the catheter's tip 1–3 cm from the wound and press the trigger. Once finished, release the trigger, wait a few seconds, then withdraw the catheter and switch off the spraying device.

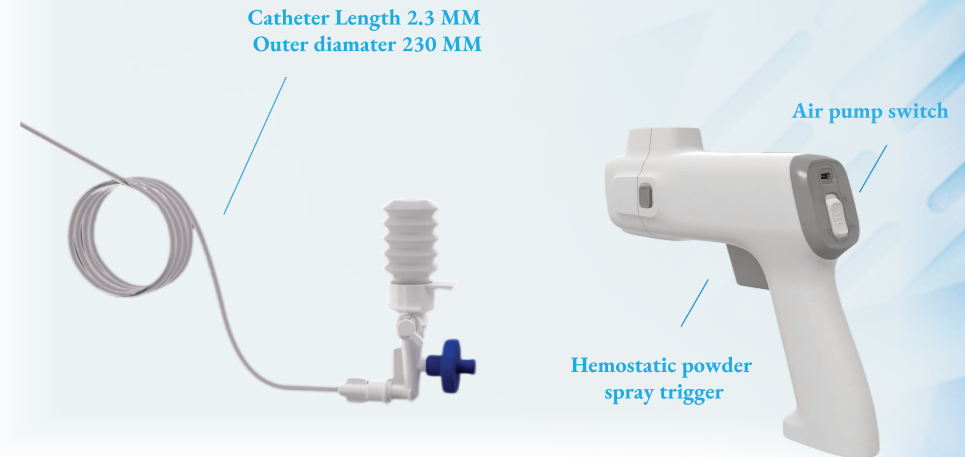

A New Choice For Hemostasis In Endoscopic Surgery

- Specially designed for Endoscopic Surgery.
- Precise powder delivery to ensures quick bleeding control
- User-friendly workflow, can be managed by one person.
- Versatile surgical unity, suitable for a wide range of other surgical procedures
- High-quality build, with rechargeable/long battery life, designs for repeated use.